

CELLUShapeTM

CAVITATION VACUUM RF

Body contouring and fat reduction



Sincoheren
www.sincoherengroup.com

PRODUCT DETAILS

The first device combining **five technologies** in one machine for the optimal **body contouring**, **skin tightening** and **weight loss** treatments, to maximize your ROI!

Intuitive Interface

The system's user interface helps optimize treatments by offering both continuous and pulse modes.



Technology : Five different technologies including IR(infrared),RF(Bi-polar radio frequency),Vacuum,Cavitation and Auto Roller Massage

Multi-treatments: Five technologies in one device to cover all face and body treatments


Handpiece : Five different treatment heads meet different functional and positional treatment needs

Handle cord: 2.0 meter treatment cord to

WORKING THEORY


CelluShape

is the 4th edition of Kuma Shape series

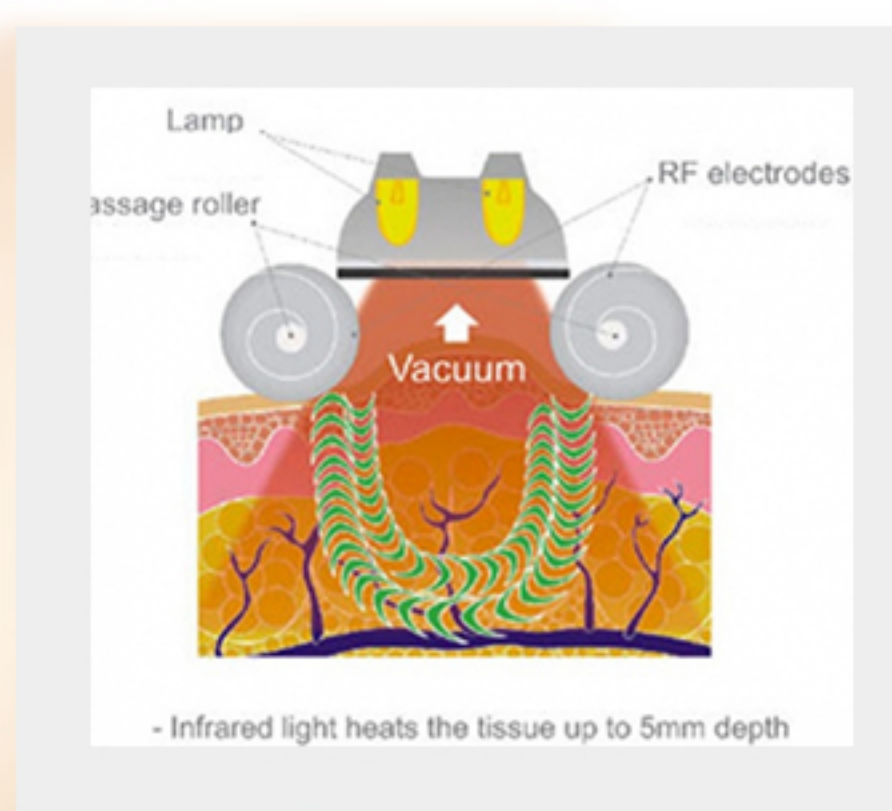

A new and promising technology for non-surgical body contouring, fat and cellulite reduction, face neck lifting and eyes wrinkle removal. It is safe, effective with worldwide proved clinical efficacy.

- Infrared light(IR) heats the tissue up to 3 mm depth.
- Bi-polar radio frequency (RF) heats tissue up to ~15mm depth.
- Vacuum+/- massage mechanisms enable precise targeting of energy to the tissue.
- 40KHZ Cavitation target and destroy subcutaneous fat cells in the fat cell membrane.

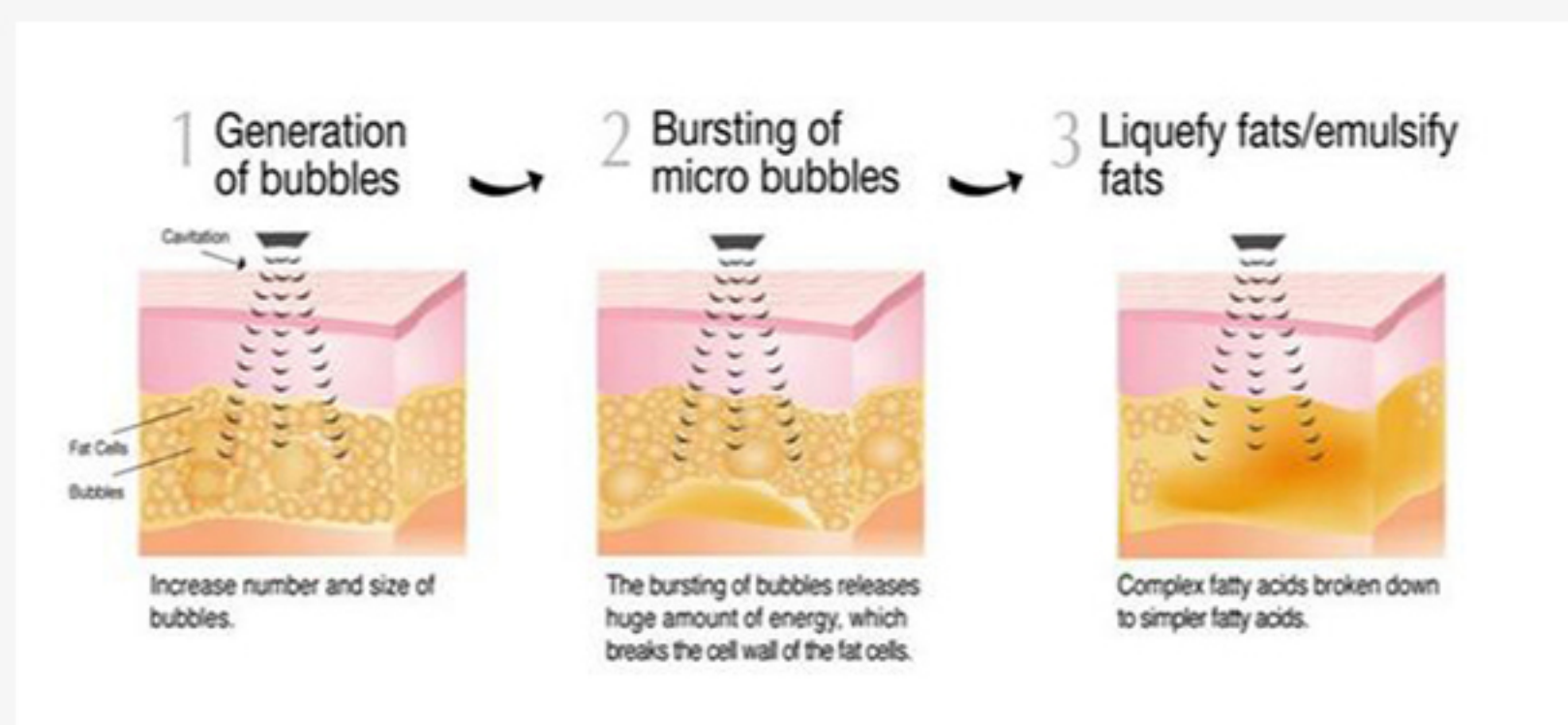

Advantages

Foldable screen :

The foldable screen design can meet the needs of treatment under different conditions

1



Machine Holder and roller design :

Humanized holder design and roller design, easy to move the machine during treatment

2



Built-in filters on the side for easy cleaning and save ongoing cost

3



APPLICATION



Face



Eyes



Back



Waist



Abdomen



Hips



Thigh



Lower leg



BEFORE & AFTER

Before


After 6 weeks



Photos courtesy of Dr. Laura Johson, France

Before

After 9 weeks

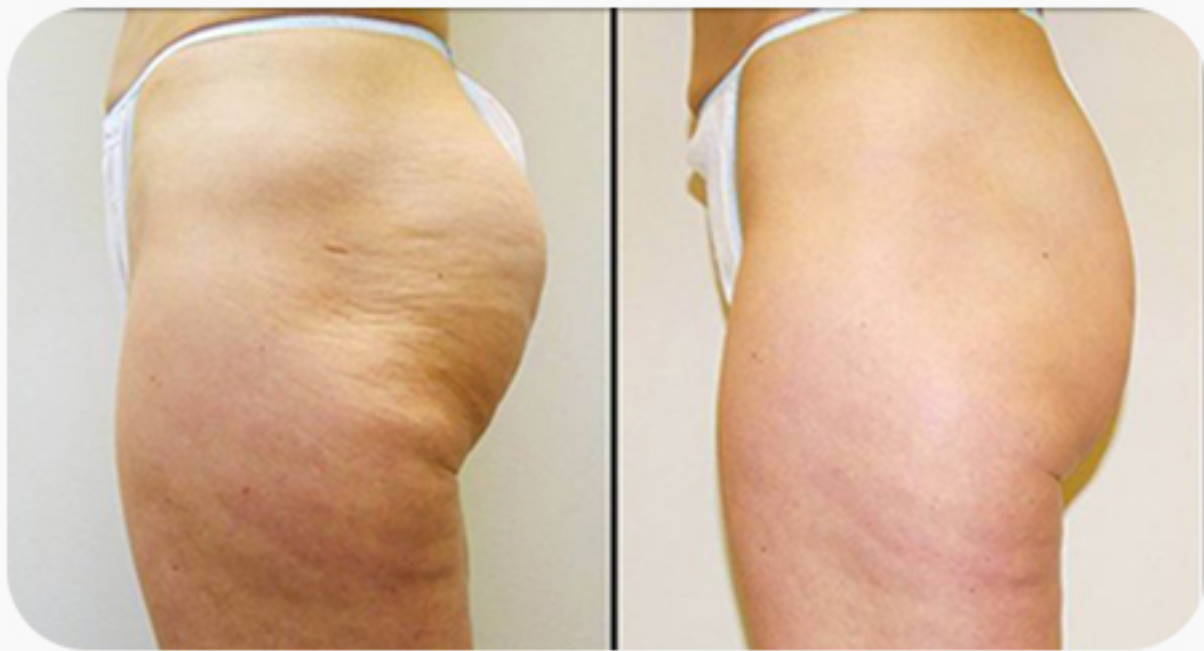

Photos courtesy of Dr. Kim Nuygen, USA

Before

After 4 weeks



Photos courtesy of Dr. Silke Mark, Holland

Before

After 8 weeks



Photos courtesy of Dr. Natasha Tang, Poland

SPECIFICATION

Items	Specification
Pulse width	0.5s-7.5s
Negative pressure	10KPA-80KPA
Rev of roller	0-36 rpm
Working mode for roller	4 types
Frequency	1-10MHz
Laser power	MAX 20W
Cavitation	40KHz
Treatment area	4mmx7mm;8mmx25mm; 30mmx50mm;40mmx60mm
Rated input power	750VA
Mode of power supply	AC230V±10%,50Hz±1Hz/ AC110V±10%,60Hz±1Hz(Optional)
G.W.	62.8Kg,16.5Kg
Packing size	70*170*62cm,67*57*31cm



Add: A-4,Sinotrans Plaza,No.43 Xizhimen North
St,Beijing 100044 C
Tel: +86 18701341926
Web: www.sincoherengroup.com

