



Sincoheren
www.sincoherengroup.com



coolplas®

Cryolipolysis


Painless and non-invasive

Exclusively treat obesity

State-of-the-art

Product details

The main interface and the handle interface display real-time working temperature, you can pay attention to the machine temperature and adjust the temperature at any time

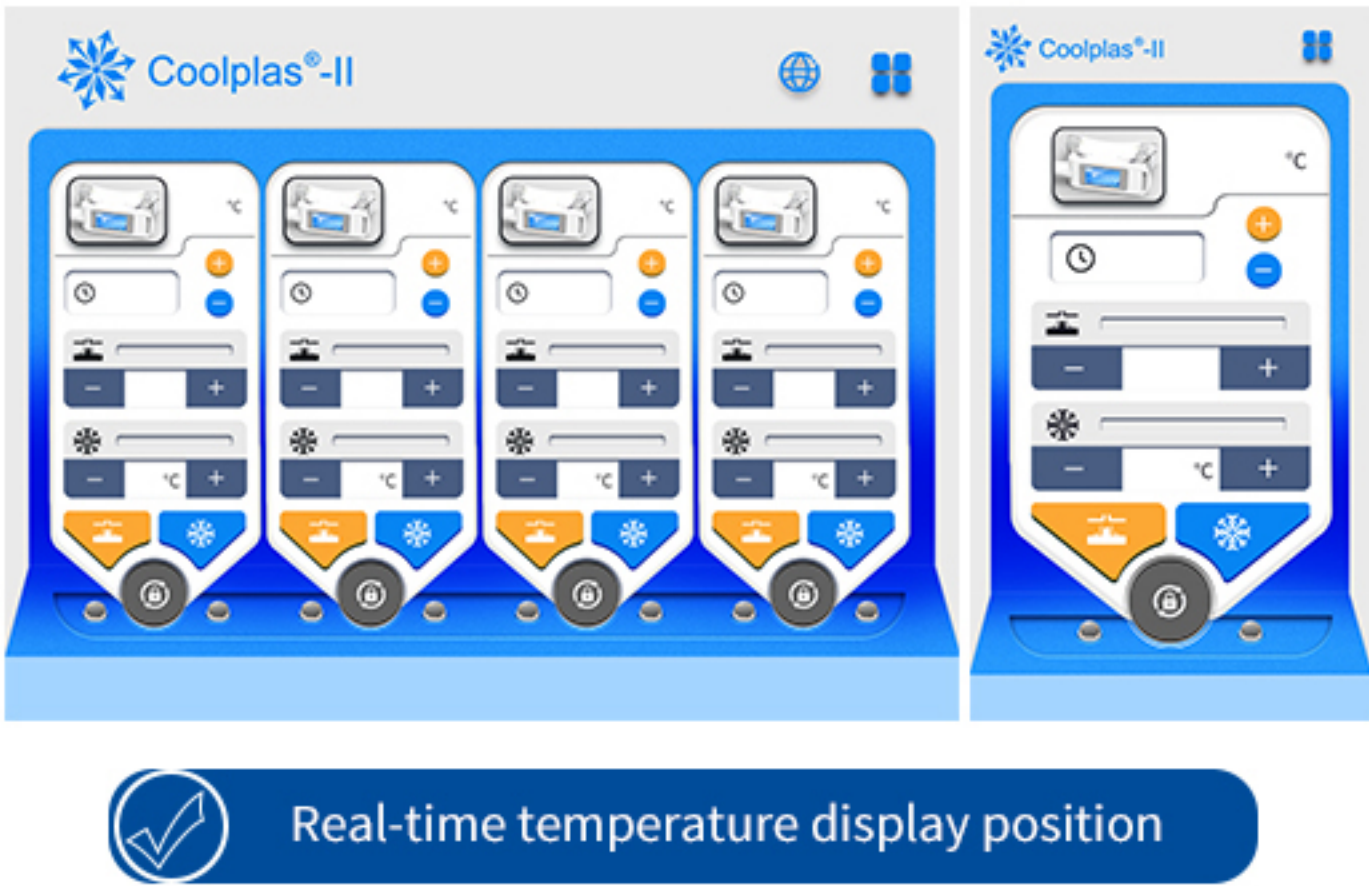


Four handles can work simultaneously and asynchronously






ICE 1.0
Cavity size:78mm*33mm
Covered area:23cm²

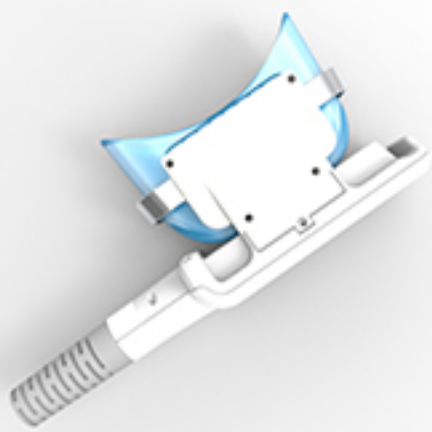


The curved cup is designed for better placement and suitable for curved or narrow body parts

ICE 3.0
Cavity size:180mm*60mm
Covered area:98cm²

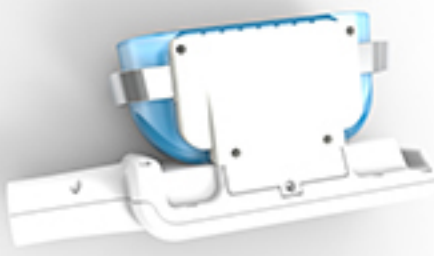


Suitable for longer vertical fat areas (ie abdomen and love handles), also recommended for arms and inner thighs

ICE 3.0F
Cavity size:170mm*60mm
Covered area:86cm²

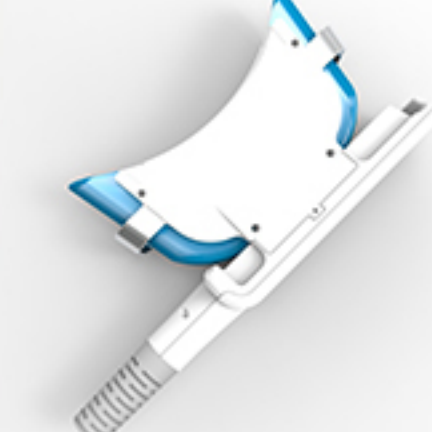

Best for belly contour. The most commonly used treatment probe in this area

ICE 6.0
Cavity size:270mm*60mm
Covered area:175cm²

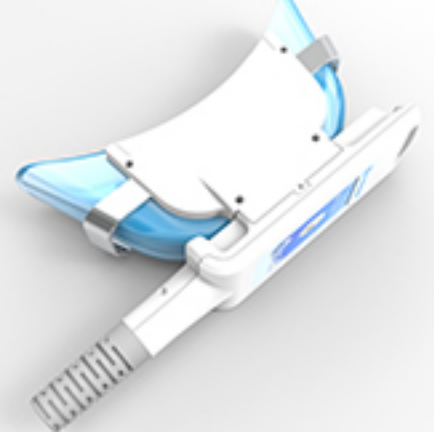

2.5 times larger than ICE3.0 treatment head. It is recommended to reduce or reduce swelling

ICE 8.0(optional)
Cavity size:340mm*60mm
Covered area:243cm²

Principle

As triglyceride in fats will be converted into solid in particular low temperatures, it uses advanced cooling technology to selectively target fat bulges and eliminate fat cells through a gradual process that does not harm the surrounding tissues, reduce unwanted fat, When fat cells are exposed to precise cooling, they trigger a process of natural removal that gradually reduces the thickness of the fat layer. And the fat cells in the treated area are gently eliminated through the body's normal metabolism process, to eliminate unwanted fat. Precisely controlled cooling is applied through the skin to the fat layer, and eliminates local fat. The whole time of procedure is about 1~2 hour



1

Many of us have bulges of stubborn fat like muffin tops

A realistic illustration of a person's midsection showing a prominent 'muffin top' bulge.

2

Those unwanted buldges contain fat cells,which can be resistant to diet and exercise but not to Kryolipolyse

A cross-sectional diagram of skin showing a layer of yellow, spherical fat cells.

3

Kryolipolyse uses controlled cooling to target and crystallize fat cells

A cross-sectional diagram of skin showing fat cells that are partially covered by a blue, crystalline substance, representing the cooling process.

4

Crystallized fat cells gradually die off,thenare naturally eliminated from your body

A cross-sectional diagram of skin showing the blue crystalline substance being removed from the area, leaving the yellow fat cells exposed.

5

In the weeks and months follwing treatment,remaining fat cells condense,reducing the fat layer

A cross-sectional diagram of skin showing a single layer of yellow fat cells, indicating a reduction in the fat layer's thickness.


6

Kryolipolyse lets you say goodbye to stubborn fat

A realistic illustration of a person's midsection showing a significantly reduced 'muffin top' bulge, indicating a slimmer waistline.

Advantage

Machine appearance structure

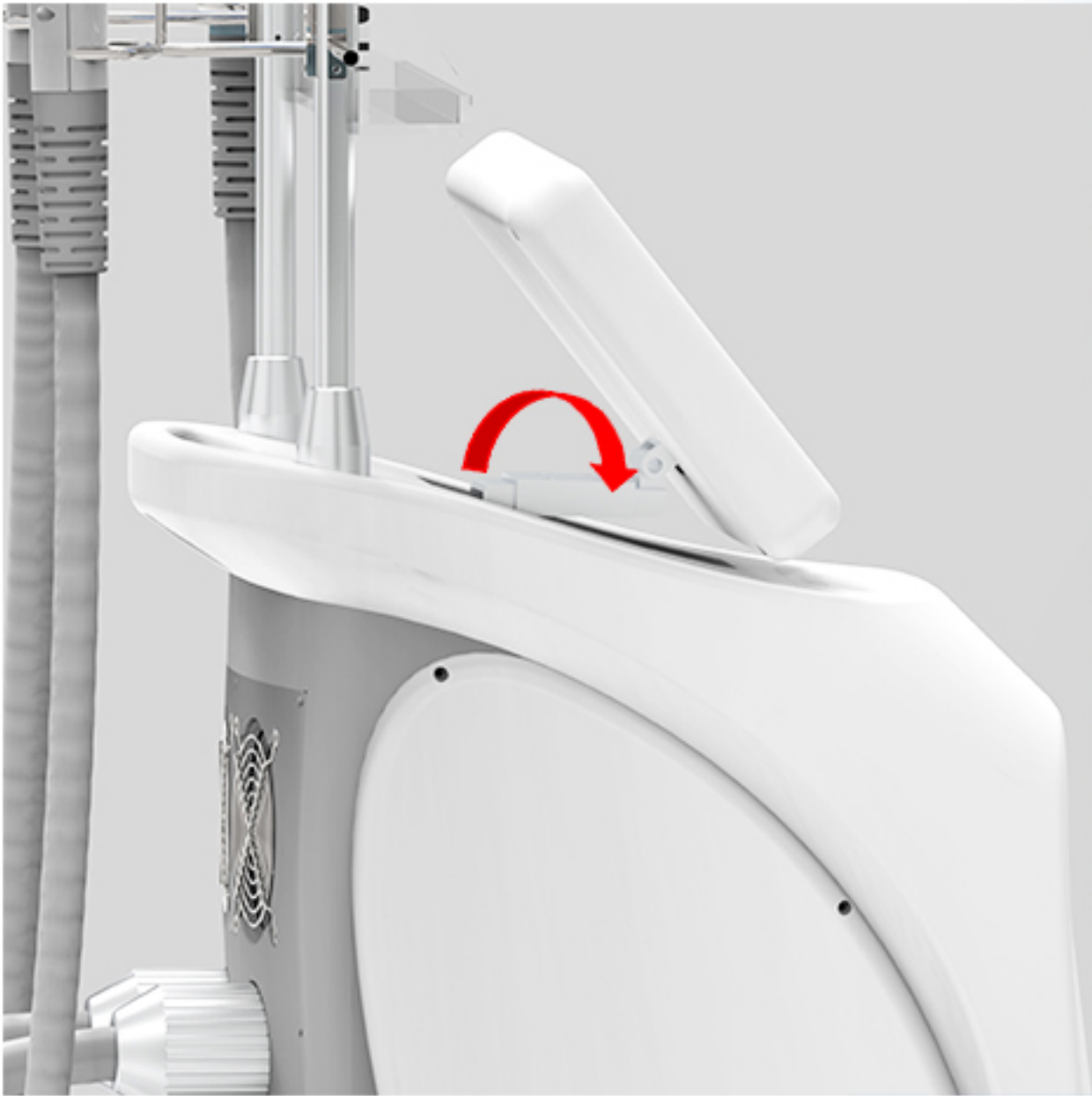

Foldable screen

The machine screen can be rotated up and down about 90 degrees, and the height and angle can be adjusted freely to adapt to different light




The machine can be opened on three sides, because the filter needs to be checked frequently, the design of the maintenance window position is reserved in advance



Filter reserved window

Application

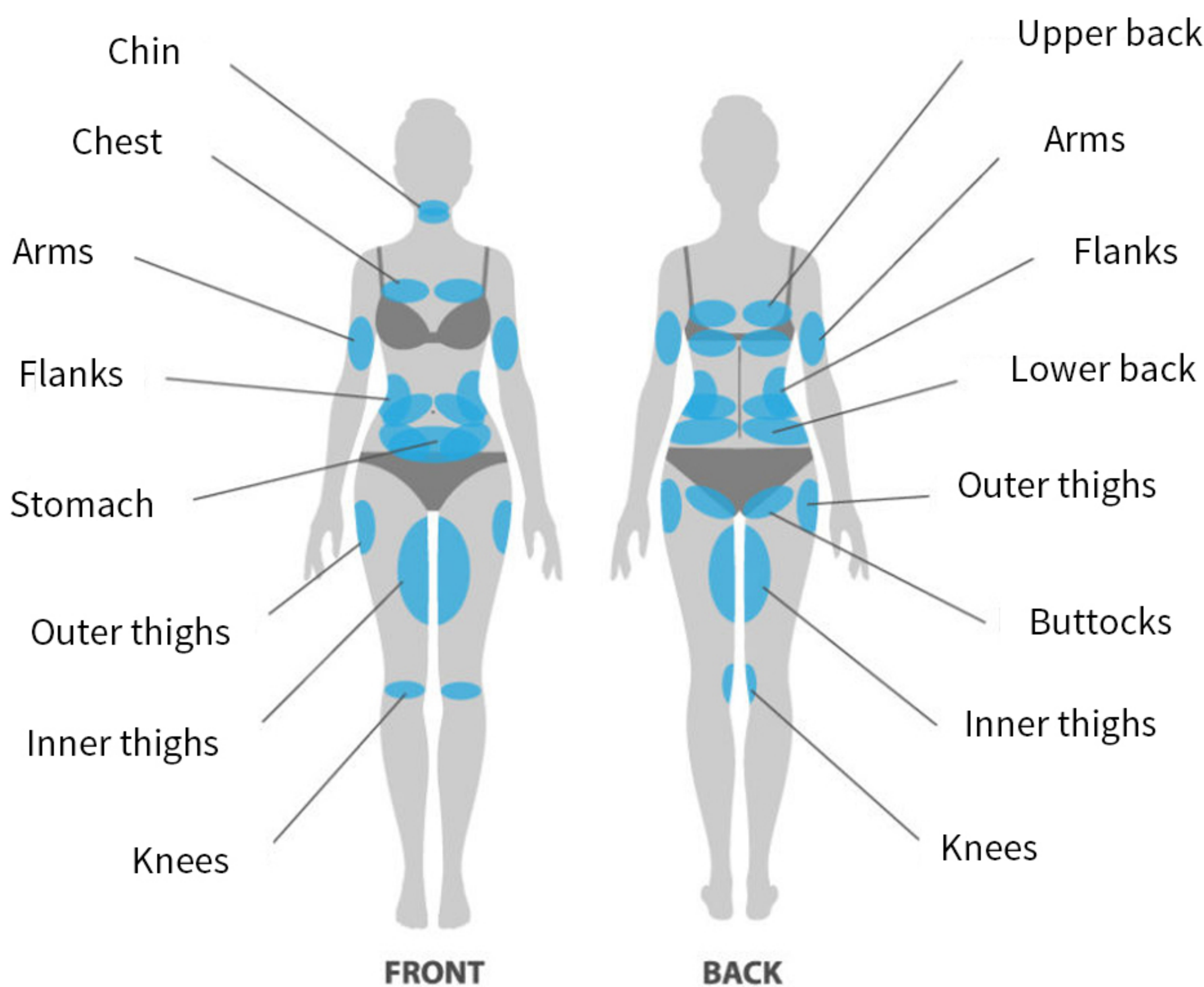

**FREEZE
STUBBORN
FAT**
no surgery, just
Coolplas


Before & After



Before

8 weeks after
First coolplas' session

9 weeks after
Second coolplas' session
Photos courtesy of Edward Becker,MD

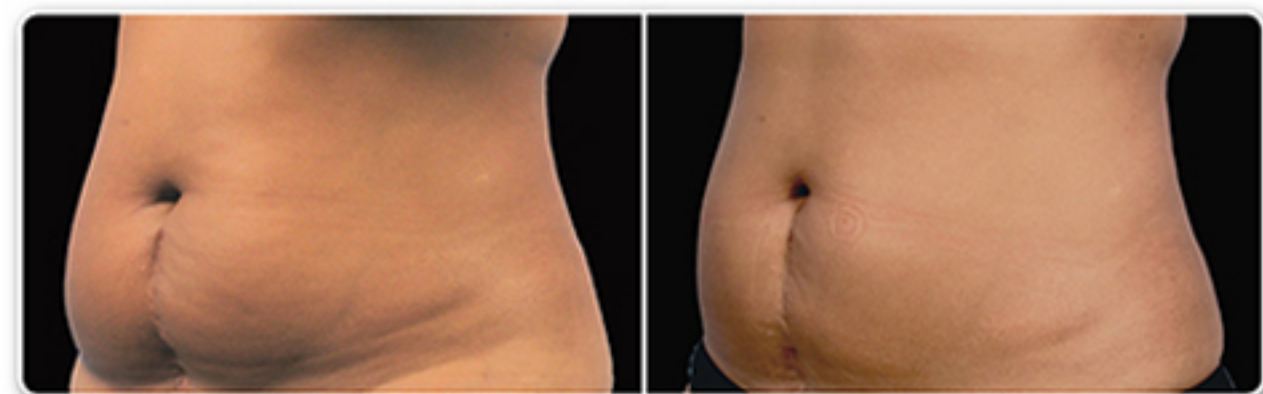

Before

60 days after
One coolplas' treatment
Photos courtesy of Charles Mok,DO



Before

60 days after
Photos courtesy of Grant Stevens,MD

Specification

Temperature rage	-8℃~ 0℃
Temperature accuracy	±0.1℃
Treatment area size	Mid:250mmx50mm Min:145mmx50mm Fit:140mmx50mm
Negative pressure	0-0.07 M Pa
Duration	1~60min
Electrical requirement	230VAC 50Hz 1500W
Dimension(WxDxH)	300mmx400mmx900mm
Weight	30Kg



Add: A-4,Sinotrans Plaza,No.43 Xizhimen North
St,Beijing 100044 C
Tel: +86 18701341926
Web: www.sincoherengroup.com

